


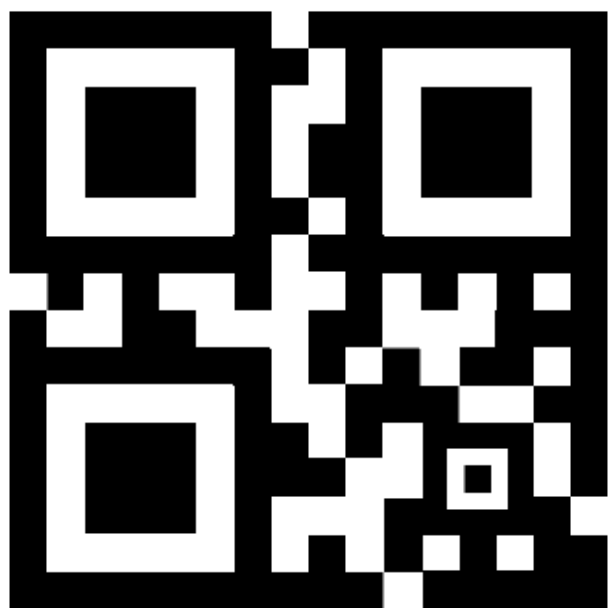

LabTrackQR

Scanner Configuration Manual

This manual contains the essential setup sequence to configure **Zenwire W213** scanners for seamless communication with **LabTrackQR**.

Instructions:

Please scan the barcodes on the following pages in the exact order they are presented. Do not skip any steps.

Important:

Please pause for a few seconds between each scan to allow the device to process and apply the configuration.


After each successful scan, the scanner will emit a beep to confirm that the setting has been applied.

Initialization & Factory Reset

To ensure seamless integration with the system, we must first clear any existing configurations and start with a clean slate.

Step 1: Enter Setup Mode


Scanning this barcode unlocks the device's configuration settings. This is always required before making any changes.

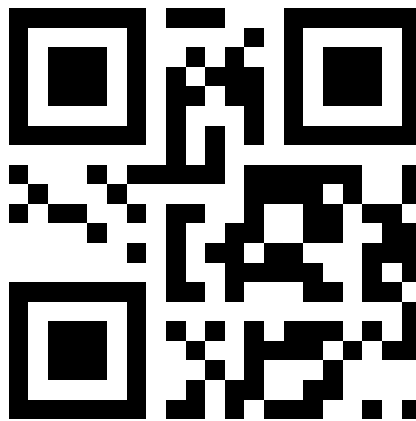

[QR CODE: S_CMD_0001]

Step 2: Factory Reset

This restores the scanner to its original factory defaults, securely erasing all previous custom settings.

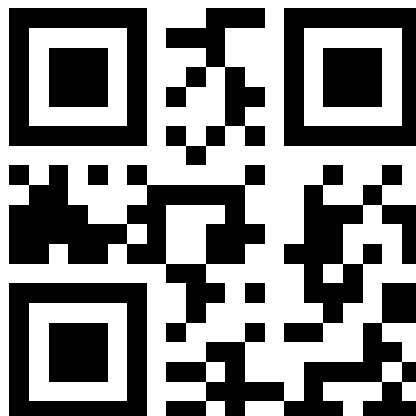


[QR CODE: S_CMD_FFFF]

Communication Protocol

LabTrackQR requires a specific communication method to correctly capture and process your scans. This section configures the scanner to use a Virtual COM Port.

Step 3: Re-Enter Setup Mode


Since the factory reset cleared the device memory, you must unlock the configuration settings once again before applying new rules.

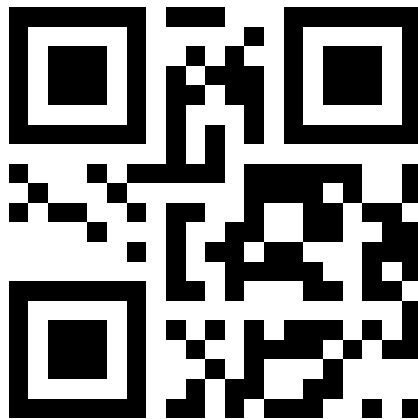

[QR CODE: S_CMD_0001]

Step 4: Enable USB Virtual Serial Port

This establishes the correct data transfer protocol between your scanner and the LabTrackQR and prevents the scanner from acting like a keyboard and typing data into active text fields.

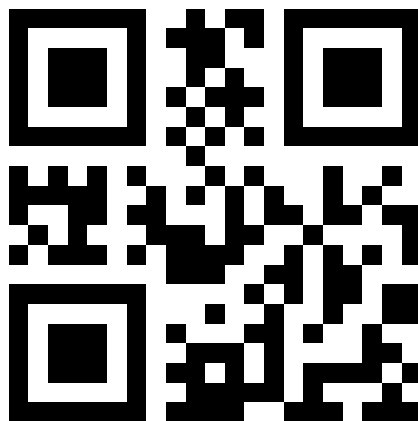


[QR CODE: S_CMD_01H3]

Scanning Mode

By default, it's highly recommend using the manual trigger mode for maximum control, precision, and to prevent accidental scans.


Step 5: Manual Mode

This setting ensures the scanner will only activate and read a barcode when you physically pull the trigger.

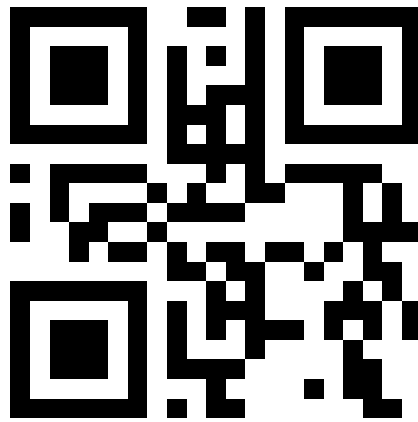

[QR CODE: S_CMD_MT00]

CONFIGURATION COMPLETE

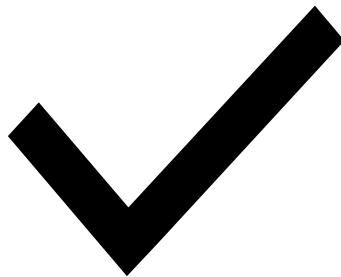

Congratulations! Your Zenwire W213 scanner is now fully optimized and ready to use.

You can now close this manual and begin scanning your items in the LabTrackQR application.


OPTIONAL: Advanced Scanning Modes

If your specific workflow requires hands-free operation, you can choose one of the alternative modes below.

Important: *If you are configuring these options much later, remember to unlock the device first by scanning QR code from **Step 1** (Page 3).*

Option A: Continuous Mode


The laser remains on constantly, allowing for rapid, continuous scanning without needing to use the trigger.

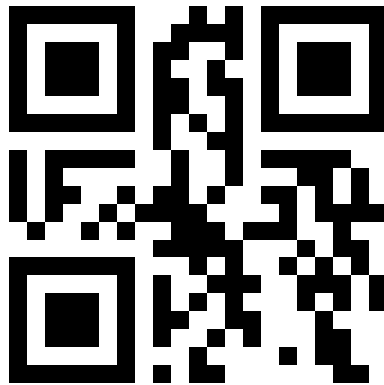

[QR CODE: S_CMD_020E]

Option B: Sensing Mode

The scanner automatically activates the laser only when it detects motion or an item placed directly under the lens and scans automatically without needing to use the trigger.

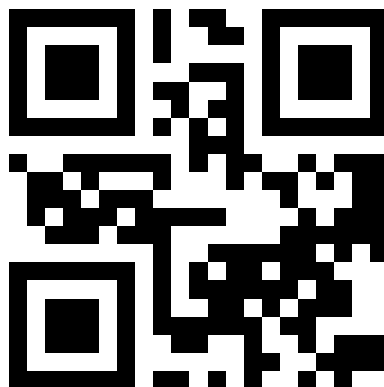

[QR CODE: S_CMD_020F]